

HARVARD UNIVERSITY
SCIENCE CENTER
PHYSICS 191/247 LABS
DATE RECEIVED: JUNE 2002
INITIALS:

Model 671 s/n 556 rev. K
Spectroscopy Amplifier
Operating and Service Manual

CONTENTS

STANDARD WARRANTY	ii
SAFETY INSTRUCTIONS AND SYMBOLS	iv
SAFETY WARNINGS AND CLEANING INSTRUCTIONS	v
1. DESCRIPTION	1
1.1. GENERAL	1
2. SPECIFICATIONS	3
2.1. PERFORMANCE	3
2.2. CONTROLS AND INDICATORS	4
2.3. INPUTS	4
2.4. OUTPUTS	5
2.5. ELECTRICAL AND MECHANICAL	5
3. INSTALLATION	6
3.1. POWER CONNECTION	6
3.2. PREAMPLIFIER CONNECTION	6
3.3. PULSED RESET PREAMPLIFIERS AND INHIBIT IN CONNECTION	6
3.4. CONNECTION OF TEST PULSE GENERATOR	6
3.5. SHAPING CONSIDERATIONS	7
3.6. LINEAR OUTPUT CONNECTIONS	7
3.7. PILE-UP REJECTION USING PUR OUTPUT	8
3.8. LIVETIME CORRECTION USING BUSY OUTPUT	8
3.9. INPUT COUNT RATE USING CRM OUTPUT	8
4. OPERATING INSTRUCTIONS	9
4.1. INITIAL TESTING AND OBSERVATION OF PULSE WAVEFORMS	9
4.2. STANDARD SETUP PROCEDURES	9
4.3. POLE-ZERO ADJUSTMENTS FOR RESISTIVE-FEEDBACK PREAMPLIFIER	10
4.4. BASELINE RESTORER (BLR) SETTING	11
4.5. INTERNAL CONTROLS	12
4.6. DIFFERENTIAL INPUT MODE	13
4.7. SYSTEM THROUGHPUT	14
4.8. CHARGE COLLECTION OR BALLISTIC DEFICIT EFFECTS	15
4.9. PILE-UP REJECTOR (PUR) AND LIVETIME CORRECTOR	16
4.10. OPERATION WITH SEMICONDUCTOR DETECTORS	17
4.11. OPERATION IN SPECTROSCOPY SYSTEMS	20
4.12. OTHER EXPERIMENTS	21
5. MAINTENANCE	25
5.1. TEST EQUIPMENT REQUIRED	25
5.2. PULSER TEST	25
5.3. SUGGESTIONS FOR TROUBLESHOOTING	26
5.4. FACTORY REPAIR	26
5.5. TABULATED TEST POINT VOLTAGES	26

ORTEC MODEL 671 SPECTROSCOPY AMPLIFIER

1. DESCRIPTION

1.1. GENERAL

The ORTEC Model 671 high-performance, energy spectroscopy amplifier is ideally suited for use with germanium, silicon surface-barrier, and Si(Li) detectors. It can also be used with scintillation detectors and proportional counters. The Model 671 input accepts either positive or negative polarity signals from a detector preamplifier and provides a positive 0 to 10-V output signal suitable for use with single- or multichannel pulse height analyzers. Its gain is continuously variable from 2.5 to 1500.

Automation of critical adjustments makes the 671 easy to set up with any detector, while minimizing the required operator expertise.

A front-panel switch on the Model 671 provides the choice of either a triangular or a Gaussian pulse shape on the UNIPOLAR output connector. (Fig. 1.1) The noise performance of the triangular pulse shape is equivalent to a Gaussian pulse shape having a 17% longer shaping time constant. In applications where the series noise component is dominant, and the pile-up rejector is utilized, the triangular shape will generally offer the same deadtime and slightly lower noise than the Gaussian pulse shape. A front-panel switch permits selection of the optimum shaping time constant for each detector and application. Six time constants in

the range of 0.5 to 10 μ s, and the TR I/GAUSS switch combine to offer 12 different shaping times. A bipolar output is also provided for measurements requiring zero cross-over timing.

To minimize spectrum distortion at medium and high counting rates (Fig. 1.2), the unipolar output incorporates a high-performance, gated, baseline restorer with several levels of automation. Automatic positive and negative noise discriminators ensure that the baseline restorer

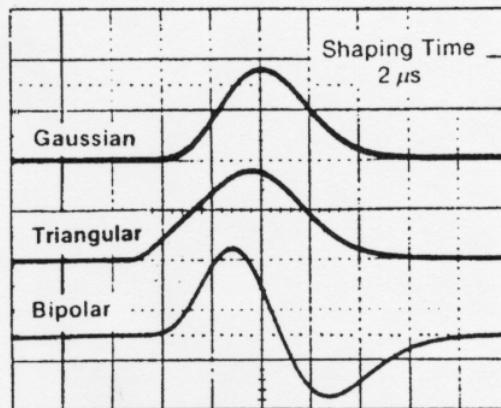


Fig. 1.1. Gaussian, Triangular, and Bipolar Pulse Shapes for a 2- μ s Shaping Time. Vertical scale, 5 V per division; horizontal scale, 2 μ s per division.

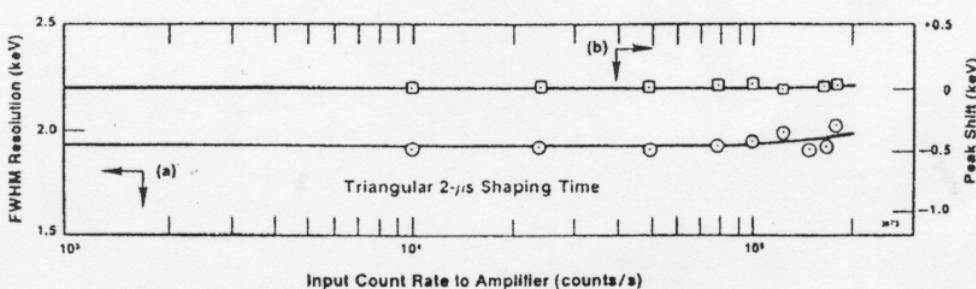


Fig. 1.2. (a) Resolution and (b) Peak Position Stability as a Function of Counting Rate.
See specifications for spectrum broadening and spectrum shift.

operates only on the true baseline between pulses in spite of changes in the noise level. No operator adjustment of the baseline restorer is needed when changes are made in the gain, the shaping time constant, or the detector characteristics. Negative overload recovery from the reset pulses generated by transistor reset preamplifiers and pulsed optical feedback preamplifiers is also handled automatically. A monitor circuit gates off the baseline restorer and provides a reject signal for a multichannel analyzer until the baseline has safely recovered from the overload.

Several operating modes are selectable for the baseline restorer. For making either a manual or automatic PZ adjustment, the PZ position is selected. This position can also be used where the slowest baseline restorer rate is desired. For situations where low frequency noise interference is a problem, the HIGH rate can be chosen. On detectors where perfect PZ cancellation is impossible, the AUTO baseline restorer rate provides the optimum performance at both low and high counting rates.

A front-panel limit (LIM) push button is included with the unipolar output to facilitate monitoring the accuracy of the PZ adjustment on an oscilloscope. When pressed, this button inserts a diode limiter in series with the unipolar output connector. This prevents overload distortions in the oscilloscope when using the more sensitive amplitude scales required for observing the PZ adjustment.

An efficient pile-up rejector is incorporated in the 671 Spectroscopy Amplifier. It provides an output logic pulse for the associated multichannel analyzer to suppress the spectral distortion caused by pulses piling up on each other at high counting rates (Fig. 1.3). The fast amplifier in the pile-up rejector includes a gated baseline restorer with its own

automatic noise discriminator. A multicolor pile-up rejector LED on the front panel indicates the throughput efficiency of the amplifier. At low counting rates the LED flashes green. The LED turns yellow at moderate counting rates and red when pulse pile-up losses are >70%.

When long connecting cables are used between the detector preamplifier output and the amplifier input, noise induced in the cable by the environment can be a problem. The Model 671 provides two solutions. For low to moderate interference frequencies the differential input mode can be used with paired cables from the preamplifier to suppress the induced noise. At high frequencies a common mode rejection transformer built into the 671 input reduces noise pick-up. The transformer is particularly effective in eliminating interference from the display raster generators in personal computers.

All toggle switches on the front panel lock to prevent accidental changes in the desired settings.

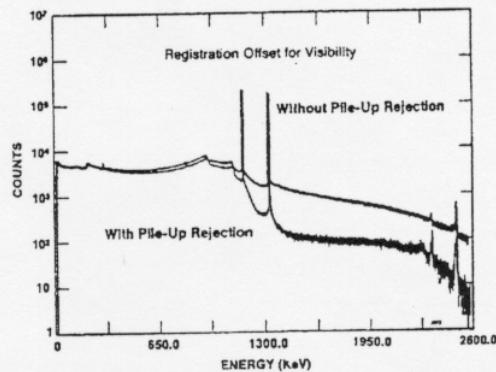


Fig. 1.3. Demonstration of the Effectiveness of the Pile-Up Rejector in Suppressing the Pile-up Spectrum. See Pulse Pile-Up Rejector specifications.

2. SPECIFICATIONS¹

2.1. PERFORMANCE

Note: Unless otherwise stated, performance specifications are measured on the unipolar output with 2- μ s Gaussian shaping, the manual PZ mode, and the AUTO BLR mode.

GAIN RANGE Continuously adjustable from 2.5 to 1500. Gain is the product of the COARSE and FINE GAIN controls.

UNIPOLAR PULSE SHAPES Switch selection of a nearly triangular pulse shape or a nearly Gaussian pulse shape at the UNI output (Fig. 1.1, Table 2.1).

BIPOLAR OUTPUT PULSE SHAPE Rise of the bipolar output pulse from 0.1% to maximum amplitude is 1.65 times selected SHAPING TIME. Zero cross-over of the bipolar output pulse delayed from the maximum amplitude of Gaussian UNIPOLAR output by 0.33 times selected SHAPING TIME.

INTEGRAL NONLINEARITY (UNIPOLAR Output) $\pm 0.025\%$ from 0 to +10 V.

NOISE Equivalent input noise <5.0 μ V rms for gains >100, and <4.5 μ V rms for gains >1000.

TEMPERATURE COEFFICIENT (0 to 50 °C)

Table 2.1. Unipolar Pulse Shape Parameters for the Triangular and Gaussian Pulse Shapes

Time Interval	Shaping Time Multiplier ^a	Triangular	Gaussian
From start of input pulse to maximum amplitude of unipolar output pulse		2.6	2.8
Rise of output pulse from 0.1% to maximum amplitude		2.4	2.0
Width of output pulse at 50% of maximum amplitude		2.5	2.0
Width of output pulse at 1% of maximum amplitude		5.6	5.0
Width of output pulse at 0.1% of maximum amplitude		6.9	6.3

^aTime interval equals the selected front-panel SHAPING TIME multiplied by the Shaping Time Multiplier.

¹ Specifications subject to change without notice.

[†] Results may not be reproducible if measured with a detector producing a large number of slow-risetime pulses or having quality inferior to the specified detector.

DIFFERENTIAL INPUT Differential nonlinearity <±0.012% from -9 V to +9 V. Maximum input ±10 V (dc plus signal). Common mode rejection ratio >1000.

PULSE PILE-UP REJECTOR

Threshold Automatically set just above noise level on fast amplifier signal. Independent of slow amplifier BLR threshold.

Minimum Detectable Signal Limited by detector and preamplifier noise characteristics.

Pulse Pair Resolution Typically 500 ns. Measured using the ^{60}CO 1.33-MeV gamma-ray under the following conditions: 10% efficiency germanium detector, 4-V amplitude for the 1.33-MeV gamma-ray at the unipolar output, 50,000 counts/s.

2.2. CONTROLS AND INDICATORS

FINE GAIN Front-panel, 10-turn precision potentiometer with locking, graduated dial provides continuously variable, direct reading, gain factor from 0.5 to 1.5.

COARSE GAIN Front-panel, eight-position switch selects gain factors of 5, 10, 20, 100, 200, 500, and 1000.

SHAPING TIME Six-position switch on the front panel selects shaping times of 0.5, 1, 2, 3, 6, and 10 μs for the pulse-shaping filter network.

MODE Two-position locking toggle switch on the front panel selects either GAUSS (Gaussian) or TRI (Triangular) pulse shaping for the UNI (unipolar) output.

INPUT POS/NEG Front-panel, two-position locking toggle switch accommodates either positive or negative input polarities.

NORM/DIFF Two-position slide switch mounted on the printed circuit board selects the normal (NORM) or differential (DIFF) input modes. In the NORM position, both front- and rear-panel INPUT connectors function as the same normal input for the preamplifier signal cable. In the DIFF mode, the rear-panel INPUT connector becomes a differential ground reference input, and the front-panel INPUT remains the normal input for the preamplifier signal cable. In the DIFF mode, the preamplifier signal cable is connected to the front-

panel INPUT and a cable having its center conductor connected to the preamplifier ground through an impedance matching resistor is connected to the rear-panel INPUT. The input impedance matching resistor must match the output impedance of the preamplifier.

BAL (Differential Input Gain Balance) A 20-turn potentiometer mounted on the PC board inside the module allows the gains of normal and differential reference inputs to be matched for maximum common mode noise rejection in DIFF mode.

PZ ADJUSTMENT 20-turn potentiometer on the front panel permits screwdriver adjustment of the PZ cancellation. The adjustment covers preamplifier exponential decay time constants from 40 μs to ∞ . For transistor reset preamplifiers or pulsed optical feedback preamplifiers, set the PZ adjustment fully counterclockwise.

LIMIT PUSHBUTTON Inserts a diode limiter in series with the front-panel UNI output connector. Prevents overload distortions in the oscilloscope when observing accuracy of the PZ adjustment on the more sensitive oscilloscope ranges.

BLR A front-panel, three-position, locking, toggle switch selects the baseline restorer rate. PZ position offers lowest fixed rate, for adjusting PZ cancellation. AUTO position matches the rate of the PZ position at low counting rates, but increases the restoration rate as the counting rate rises. HIGH rate position is provided for suppressing low frequency interference.

PUR ACCEPT/REJECT LED Multicolor LED indicates percentage of pulses rejected because of pulse pile-up. LED appears green for 0-40%, yellow for 40-70%, and red for >70% rejection.

2.3. INPUTS

INPUT (Front Panel) Front-panel, BNC connector accepts preamplifier signals of either polarity with risetimes less than the selected SHAPING TIME, and exponential decay time constants from 40 μs to ∞ . For the NEG INPUT switch setting, the input impedance is 1000 Ω on a coarse gain of 5, and 465 Ω at coarse gain settings ≥ 10 . For the POS INPUT switch setting, the input impedance is 2000 Ω for a coarse gain of 5, and 1460 Ω for coarse gains ≥ 10 . Input is dc-coupled, and protected to ± 25 V.

INPUT (Rear Panel) BNC connector, identical to front-panel INPUT when PWB-mounted NORM/DIFF slide switch is in the NORM position. When operating in the differential input mode with the slide switch set to DIFF, the rear-panel INPUT is used for the preamplifier ground reference connection. For the DIFF and POS INPUT switch setting, the input impedance is $1000\ \Omega$ on a coarse gain of 5, and $465\ \Omega$ at coarse gain settings ≥ 10 . For the DIFF and NEG INPUT switch setting, the input impedance is $2000\ \Omega$ for a coarse gain of 5, and $1460\ \Omega$ for coarse gains ≥ 10 . Input is dc-coupled, and protected to $\pm 25\ V$.

INH IN Rear-panel BNC input connector accepts reset signals from transistor reset preamplifiers or pulsed optical feedback preamplifiers. Positive NIM standard logic pulses or TTL levels can be used. Logic is selectable as active high or active low via printed circuit board jumpers. Inhibit input initiates the protection against distortions caused by the preamplifier reset. This includes turning off the baseline restorers, monitoring the negative overload recovery at the unipolar output, and generating PUR (reject) and BUSY signals for the duration of the overload. The PUR and BUSY logic pulses are used to prevent analysis and correct for the reset deadtime in the associated ADC or multichannel analyzer.

2.4. OUTPUTS

UNI Front- and rear-panel BNC connectors provide positive, unipolar, shaped pulses with a linear output range of 0 to $+10\ V$. Front-panel output impedance $<1\ \Omega$. Rear-panel output impedance selectable for either $<1\ \Omega$ or $93\ \Omega$ using a printed circuit board jumper. Outputs are dc-restored to $0 \pm 5\ mV$ and short-circuit protected.

Bi Front- and rear-panel BNC connectors provide bipolar shaped pulses with the positive lobe leading. The linear output range is 0 to $\pm 10\ V$. Front-panel output impedance $<1\ \Omega$. Rear-panel output impedance selectable for either $<1\ \Omega$ or $93\ \Omega$ using a printed circuit board jumper. Baseline between pulses has a dc level of $0 \pm 10\ mV$. Short-circuit protected.

CRM The Count Rate Meter output has a rear-panel BNC connector and provides a 250-ns-wide, $+5\ V$ logic signal for every linear input pulse that exceeds the pile-up inspector threshold. Output impedance is $50\ \Omega$.

BUSY Rear-panel BNC connector provides a $+5\ V$ logic pulse for the duration that the linear signals exceed the positive or negative baseline restorer thresholds, or the pile-up inspector threshold, or for the duration of the INH IN input signal. Useful for deadtime corrections with an associated ADC or multichannel analyzer. Positive NIM standard logic pulse is selectable as active high or active low via a printed circuit board jumper. Output impedance is $50\ \Omega$.

PUR Pile-Up Reject output is a rear-panel, BNC connector. Provides a $+5\ V$ NIM standard logic pulse when pulse pile-up is detected. Output also present for a pulsed reset preamplifier during reset, and reset overload recovery. Output pulse is selectable as active high or active low by means of a printed circuit board jumper. Output impedance is $50\ \Omega$. Used with an associated ADC or multichannel analyzer to prevent analysis of distorted pulses.

PREAMP Rear-panel standard ORTEC connector (Amphenol 17-10090) provides power for the associated preamplifier. Mates with power cords on all standard ORTEC preamplifiers.

2.5. ELECTRICAL AND MECHANICAL

POWER REQUIRED The Model 671 derives its power from a NIM Bin supplying $\pm 24\ V$ and $\pm 12\ V$, such as the ORTEC Model 4001A/4002A Bin/Power Supply. The power required is $+24\ V$ at $100\ mA$, $-24\ V$ at $200\ mA$, $+12\ V$ at $325\ mA$, and $-12\ V$ at $180\ mA$.

WEIGHT

Net $1.5\ kg$ (3.3 lb).

Shipping $3.1\ kg$ (7.0 lb).

DIMENSIONS Standard single-width module, $3.45\ X 22.13\ cm$ (1.35 X 8.714 in.) Front panel per DOE/ER- 0457T.

3. INSTALLATION

3.1. POWER CONNECTION

The 671 operates on power that must be provided by a NIM-standard bin and power supply such as the ORTEC 4001/4002 series. Convenient test points on the power supply control panel should be used to check that the dc voltage levels are not overloaded. The bin and power supply is designed for relay rack mounting. If the equipment is rack mounted, be sure that there is adequate ventilation to prevent any localized heating of the components that are used in the 671. The temperature of the equipment mounted in racks can easily exceed the maximum limit of 50°C unless precautions are taken.

3.2. PREAMPLIFIER CONNECTION

The Preamp connector of this amplifier is directly compatible with ORTEC preamplifiers as well as with standard Aptec, Canberra, PGT, and Tennelec (serial numbers greater than 2000) preamplifiers. Preamplifier power at +24 V, -24 V, +12 V, and -12 V is available through the Preamp connector on the rear panel.

When a BNC cable longer than ten feet is used to connect the preamplifier output to the amplifier input, the characteristic impedance of the cable should match the impedance of the preamplifier output. All ORTEC preamplifiers contain series terminations that are either 93 Ω or variable; coaxial cable type RG-62/U or RG-71/U is recommended.

3.3. PULSED RESET PREAMPLIFIERS AND INHIBIT IN CONNECTION

The 671 Amplifier is directly compatible with most pulsed reset preamplifiers such as the ORTEC TRP (Transistor Reset Preamplifier) Series. The amplifier automatically senses preamplifier resets and gates off the amplifier's baseline restorer. Preamplifier inhibit signals are not required for

proper amplifier operation; however, since the preamplifier resetting process is nonlinear by nature, spurious phantom peaks may show up in the spectra if the inhibit signal from the preamplifier is not used.

INHIBIT IN CONNECTION Connection of the PREAMPLIFIER INHIBIT OUT signal to the rear-panel INHIBIT IN connector will result in the system being disabled during the reset period and thus avoid spurious peaks in the spectra. Preamplifiers with an Inhibit time switch such as ORTEC PLUS Detector with series 132 Preamplifier can be set to position "1", which is the shortest preamp inhibit blocking time.

PZ SETTING The Amplifier's PZ control should be set fully counterclockwise (CCW) when used with a pulsed reset preamplifier.

3.4. CONNECTION OF TEST PULSE GENERATOR

THROUGH A PREAMPLIFIER The satisfactory connection of a test pulse generator such as the ORTEC 419 or 448 Pulse Generator or equivalent depends primarily on two considerations: the preamplifier must be properly connected to the 671 as discussed in Sections 3.2 and 3.3, and the proper signal simulation must be applied to the preamplifier. To ensure proper input signal simulation, refer to the instruction manual for the particular preamplifier being used.

DIRECTLY INTO THE 671 The ORTEC test pulse generators are designed for direct connection. When any one of these units is used, it should be terminated with a 100-Ω terminator at amplifier input or be used with at least one of the output attenuators set at In.

SPECIAL CONSIDERATIONS FOR POLE-ZERO CANCELLATION When a tail pulser is connected directly to the amplifier input, the Pole-Zero should be adjusted. See Section 4.3 for the pole-zero

adjustment. If a preamplifier is used and a tail pulser is connected to the preamplifier test input, it is not possible to adjust the pole-zero for both the preamplifier pole and the pole from the pulser tail.

3.5. SHAPING CONSIDERATIONS

The Shaping Time switch on the front panel of the 671 can be set to select time constants in steps of 0.5, 1, 2, 3, 6, and 10 μs . Choice of triangular and Gaussian filters doubles the time constants available for optimum resolution. Triangular shaping will usually give better results. The choice of the proper shaping time is generally a compromise between operating at a shorter time constant for accommodation of high counting rates and operating with a longer time constant for a better signal-to-noise ratio. Since the full amplitude of the preamplifier output pulse must be preserved, the peaking time (measurement time) must be large compared to preamplifier output pulse risetime. The amplifier shaping time should be greater than five times the charge collection time of the detector. Use the detector manufacturer's suggested shaping times as a starting point and adjust the shaping as your needs for resolution versus count rate vary.

GERMANIUM DETECTORS Shaping times for high-purity germanium (HPGe) detectors will vary from 1 to 6 μs using the unipolar output, depending on the size, configuration, and charge collection time of the specific detector and preamplifier. Coaxial detectors have significant variations in charge collection times due to their large volumes. Compromises must often be made since the shaping time that will give the best resolution will usually be longer than the optimum time needed for the best throughput at high counting rates.

Planar detectors require shaping times in the range of 3 to 10 μs for optimum resolution. Lithium-drifted silicon detectors, Si(Li), have similar shaping time requirements.

SILICON CHARGED PARTICLE DETECTORS These detectors have very fast risetimes on the order of 10 ns or less. A unipolar output and a 0.5- to 2- μs shaping time will generally provide optimum resolution.

SCINTILLATION DETECTORS The energy resolution of scintillation counters depends largely on the scintillator and photomultiplier, and therefore a shaping time of five times the decay-time constant of the scintillator is a reasonable choice. For NaI detectors that have a decay time constant of about 230 ns, the optimum shaping time is 1 μs . The bipolar output can be used to reduce overload effects and microphonics without sacrificing resolution.

GAS PROPORTIONAL COUNTERS Proportional counters have both short and long components in their charge collection times. The components typically fall in the 0.5- to 5- μs range, and lead to variable amounts of preamplifier output signal being lost as the amplifier shaping time constant is changed. Selection of longer shaping times ($>2 \mu\text{s}$) helps to minimize the problem caused by long risetimes. Due to the multiple components in the charge collection time, the correct pole-zero cancellation is not possible. This will often cause an undershoot if the Unipolar output is used. Bipolar shaping can be used to reduce this effect with little change in the resolution.

3.6. LINEAR OUTPUT CONNECTIONS

Since the 671 unipolar output is normally used for spectroscopy, the 671 is designed with a great amount of flexibility for the pulse to be interfaced with an analyzer. To minimize spectrum distortion at medium and high counting rates, the unipolar output incorporates a high-performance, gated baseline restorer with automatic setup. Automatic positive and negative noise discriminators ensure that the baseline restorer operates only on the true baseline between pulses in spite of changes in the noise level. For pulse-height analysis, the unipolar output must be directly connected to the input of a multichannel analyzer.

The bipolar output, with its symmetry about the baseline, can be used for cross-over timing or may be preferred for spectroscopy when operating into ac-coupled systems at high counting rates. Typical system block diagrams for a variety of experiments are described in Section 4.

3.7. PILE-UP REJECTION USING PUR OUTPUT

The PUR (Pile-Up Reject) output on the rear panel is used at the gate or pile-up reject input of a multichannel analyzer to suppress pile-up in the recorded spectrum. The fast amplifier in the pile-up rejector includes a gated baseline restorer with an automatic noise discriminator to eliminate the need for any operator adjustments. When pileup occurs, a logic true pulse is generated which lasts until the unipolar output returns to the baseline, normally a width of six times the shaping time. If used with a pulsed reset preamplifier, this output also includes a reject during the reset and recovery interval.

3.8. LIVETIME CORRECTION USING BUSY OUTPUT

The signal from the rear-panel Busy output connector provides a nominally +5 V logic pulse for the duration that the Unipolar output pulse exceeds

the baseline restorer threshold or pile-up inspector threshold or when the external INH IN is true. For livetime correction, Busy should be connected to the Busy In connector on the MCA. For optimal livetime correction with ORTEC analyzers like the ADCAM®, an internal jumper in the amplifier should be set to match the unipolar, triangular, or Gaussian mode. The output is internally jumper selectable as active low or active high. It is shipped as active high.

3.9. INPUT COUNT RATE USING CRM OUTPUT

A positive logic pulse is generated for each 671 input pulse that exceeds the pile-up inspector threshold level. The pulses are available through the CRM (Count Rate Meter) output on the rear panel and are intended for use in a count rate meter or counter to monitor the true input count rate into the amplifier.

4. OPERATING INSTRUCTIONS

4.1. INITIAL TESTING AND OBSERVATION OF PULSE WAVEFORMS

Refer to Section 6 for information on testing performance and observing waveforms using a pulser. Figure 4.1 shows some typical unipolar Gaussian, unipolar triangular, and bipolar output waveforms.

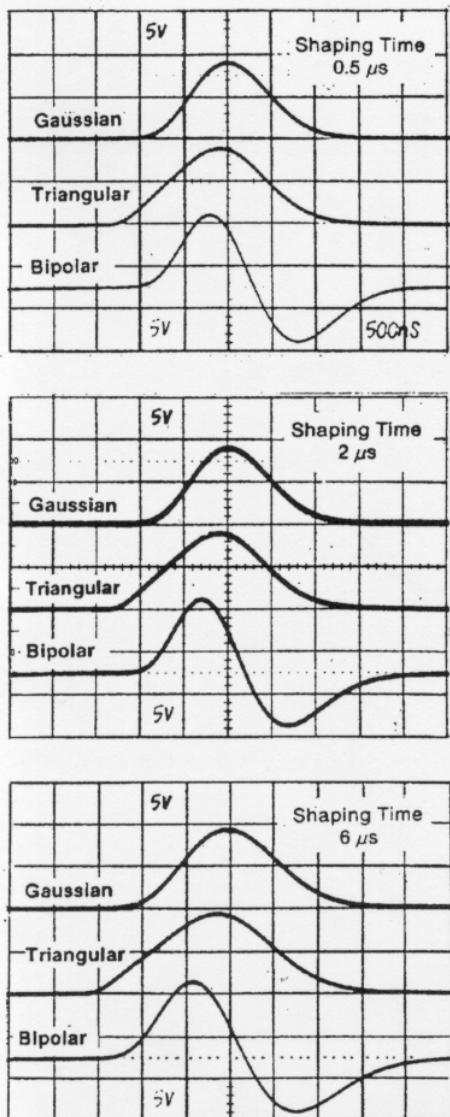


Fig. 4.1. Typical Effects of Shaping-Time Selection on Gaussian, Triangular, and Bipolar Output Waveforms.

4.2. STANDARD SETUP PROCEDURES

a. Connect the detector, preamplifier, high-voltage power supply, and amplifier into a basic system and connect the amplifier unipolar output to an oscilloscope. Connect the preamplifier power cable to the Preamp power connector on the rear panel of the 671. Turn on power in the bin and power supply and allow the electronics of the system to warm up and stabilize.

A block diagram of a typical ORTEC gamma-ray spectroscopy system is shown in Figure 4.2.

b. Set the 671 controls initially as follows:

Shaping Time	3 or 6 μ s
Mode	Triangle
Coarse Gain	20
Fine Gain	1.00
BLR	PZ
Polarity	Match preamplifier output polarity

c. Use a ^{60}Co calibration source; set about 25 cm from the active face of the detector. The unipolar output pulse from the 671 should be about 8 V, using a detector that has a preamp with a conversion gain of 300 mV/MeV.

d. Readjust the Gain control so that the higher peak from the ^{60}Co source (1.33 MeV) provides an amplifier output at about 9 V.

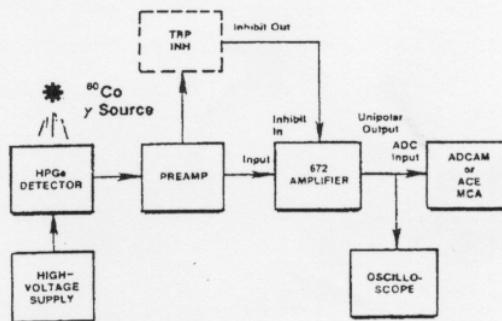


Fig. 4.2. Typical Gamma-Ray Spectroscopy System.

4.3. POLE-ZERO ADJUSTMENTS FOR RESISTIVE-FEEDBACK PREAMPLIFIER

The pole-zero adjustment is critical for good performance at high count rates in unipolar operation and for correct operation of the BLR circuit. This adjustment should be checked carefully for the best possible results. Whenever the shaping time is changed, the pole-zero must be adjusted. The bipolar output resolution is not as sensitive to misadjusted PZ, but it is important for recovery from very large overload pulses. When using a transistor reset-type preamplifier, the PZ should be set to full counterclockwise.

- Adjust the radiation source spacing from the detector to provide a count rate between 1 and 10 kHz.
- Observe the unipolar output with an oscilloscope. Increase the scope input sensitivity to 20-100 mV per vertical division. Depress the front-panel LIM push-button to limit the voltage applied to the oscilloscope. Adjust the PZ adjust control so that the trailing edge of the pulses returns to the baseline without overshoot or undershoot (Fig. 4.3). A slight bias toward an undershoot often gives the best results.

The oscilloscope used must be dc-coupled and must not contribute distortion in the observed waveforms. Oscilloscopes such as Tektronix models 465, 475, and 7904 will overload for a 10-V signal when the vertical sensitivity is < 100 mV/Div. The LIM push-button switch inserts a diode limiter in series with the front-panel UNI output connector to prevent overloading the input of the oscilloscope.

USING SQUARE WAVE THROUGH PREAMPLIFIER TEST INPUT

A more precise pole-zero adjustment of the amplifier can be obtained by using a square wave signal as the input to the preamplifier. Many oscilloscopes include a calibration output on the front panel, and this is a good source of square wave signals at a frequency of about 1 kHz. The amplifier differentiates the signal from the preamplifier so that it generates output signals of

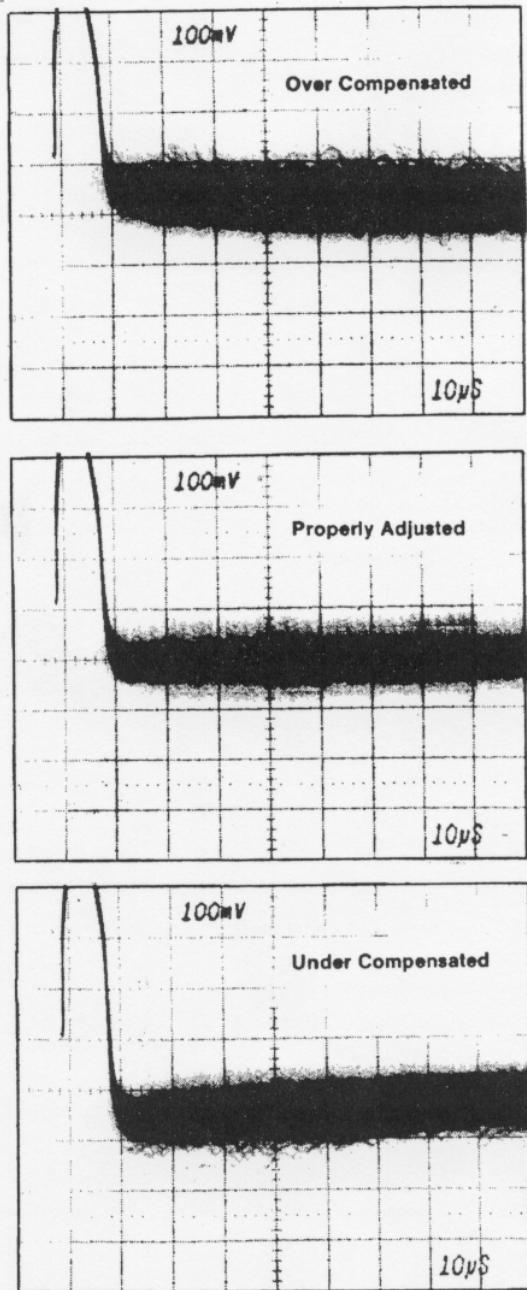


Fig. 4.3. Typical Waveforms Illustrating Pole-Zero Adjustment Effects; Oscilloscope Trigger, Busy Output; ^{60}Co Source with 1.33-MeV Peak Adjusted -9 V; Count Rate, 3 kHz; Shaping Time Constant, 2 μs .

alternate polarities on the leading and trailing edges of the square wave input signal, and these can be compared as shown in Fig. 4.4 to achieve excellent pole-zero cancellation.

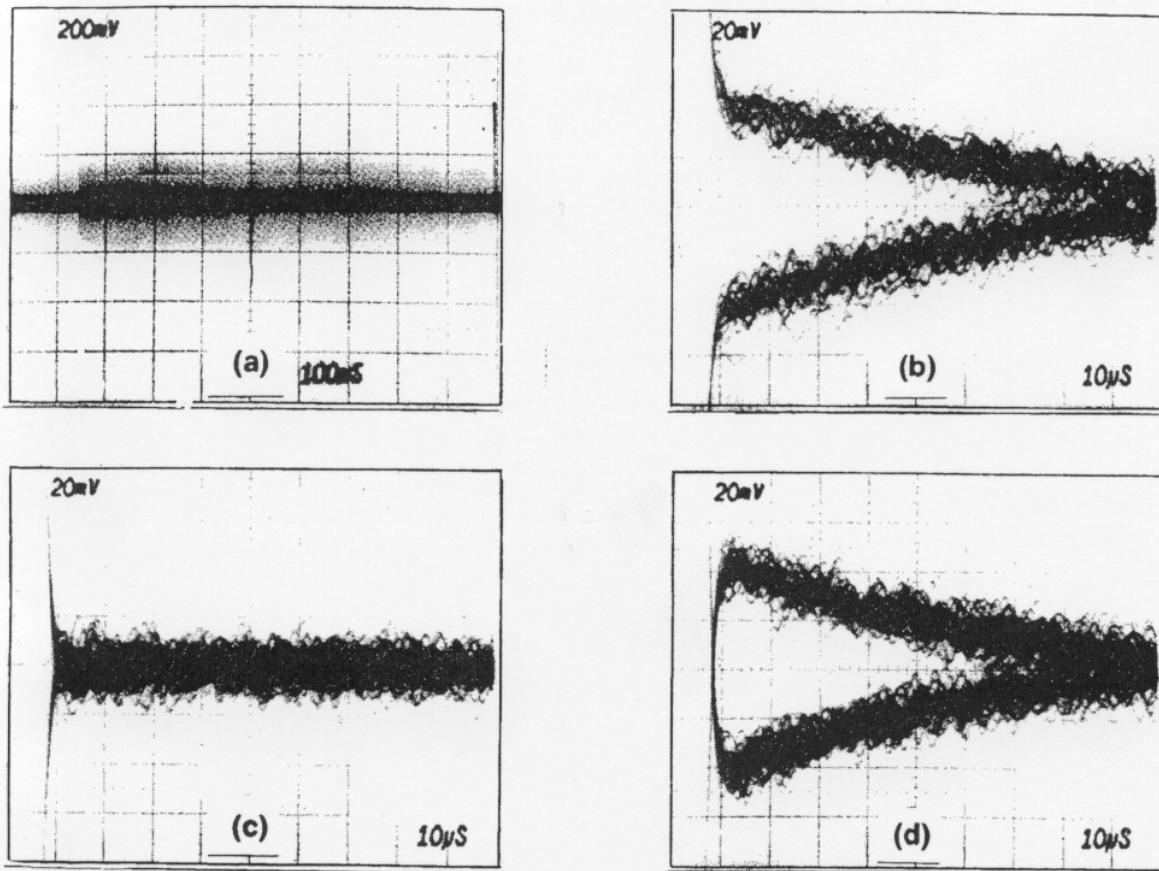


Fig. 4.4. Pole-Zero Adjustment Using a Square Wave Input to the Preamplifier. (a) PZ properly adjusted; slow trigger to separate pulses. (b) Overcompensated; fast trigger to superimpose pulses. (c) Properly adjusted; pulses superimposed. (d) Undercompensated; pulses superimposed.

Use the following procedure:

- Remove all radioactive sources from the vicinity of the detector. Set up the system as for normal operation, including detector bias.
- Set the amplifier controls as for normal operations; this includes gain, shaping, and input polarity.
- Connect the source of 1 kHz square waves through an attenuator to the Test input of the preamplifier. Adjust the attenuator so that the amplifier output amplitude is 8 to 10 volts.

d. Observe the unipolar output of the amplifier with an oscilloscope triggered from the amplifier Busy output. Adjust the PZ control for proper response according to Fig. 4.4. Depress the LIMIT push-button on the 671 while observing the adjustment on the oscilloscope display.

Figure 4.4.(a) shows the amplifier output as a series of alternate positive and negative shaped pulses. In Fig. 4.4.(b)-(c), the oscilloscope was triggered to show both positive and negative pulses simultaneously. These pictures show more detail to aid in proper adjustment.

4.4. BASELINE RESTORER (BLR) SETTING

To minimize spectrum distortion at medium and high counting rates, the unipolar output incorporates a high-performance, gated, baseline restorer with several levels of automation. Automatic positive and negative noise discriminators ensure that the baseline restorer operates only on the true baseline between pulses in spite of changes in the noise level. No operator adjustment of the baseline restorer is needed when changes are made in the gain, the shaping time constant, or the detector characteristics. Negative overload recovery from the reset pulses generated by transistor reset preamplifiers and pulsed optical feedback preamplifiers is also handled automatically to eliminate the need for operator adjustments. A monitor circuit gates off the baseline restorer and provides a reject signal for a multichannel analyzer until the baseline has safely recovered from the overload.

BLR RATE For making pole-zero adjustments, the PZ position is selected. This position can also be used where the slowest baseline restorer rate is desired.

With the BLR Rate set to AUTO, the BLR is automatically set for optimum performance throughout the usable input range for the shaping selected.

The HIGH rate can be used for situations where low or medium frequency noise interference is present and is independent of the counting rate. The HIGH rate setting is normally not used since there will be a small loss of resolution due to increased noise when used in high resolution systems.

4.5. INTERNAL CONTROLS

These controls are on the printed wiring board (PWB) and can be accessed by removing the right side cover. Figure 4.5 shows the location of these controls.

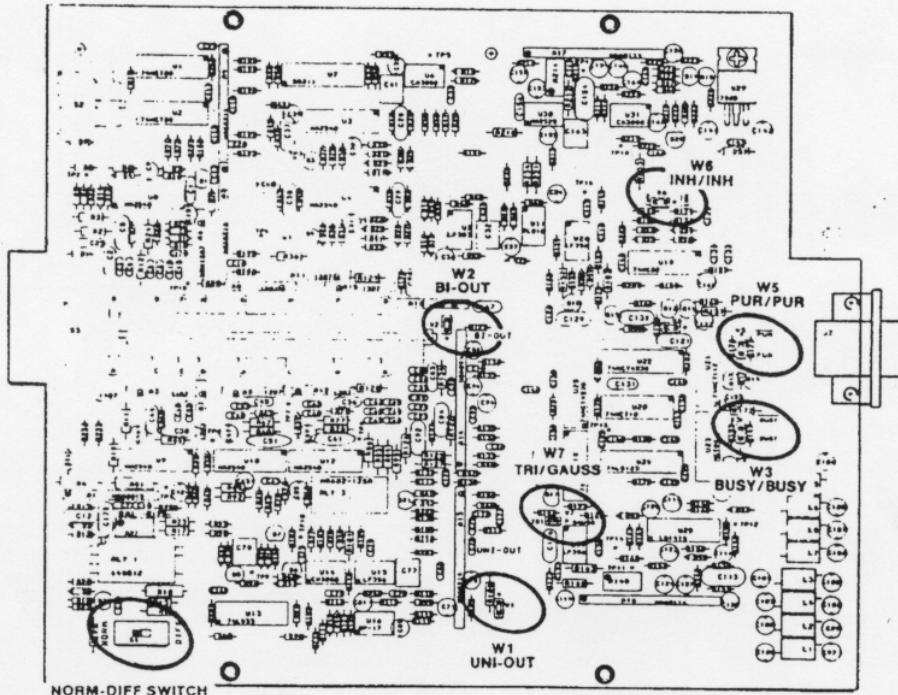


Fig. 4.5. Position of Internal Controls.

NORM-DIFF Internal PWB mounted, two-position slide switch. NORM position selects single ended inputs from front-panel input or rear-panel input connectors. In the DIFF position, the front-panel input is connected to the preamplifier signal cable, and a cable connected to the preamplifier ground through an impedance matching resistor is connected to the rear-panel input.

BAL (DIFFERENTIAL INPUT GAIN BALANCE) Internal PWB 20-turn screwdriver potentiometer allows maximization of noise rejection when using differential input. See Section 4.6.

UNI-OUT (UNIPOLAR Z_{out}) Jumper plug, W1, provides $Z_{out} \leq 1 \Omega$ or $\sim 93 \Omega$ for the rear-panel Unipolar output. Shipped in the $93-\Omega$ position.

BI-OUT (BIPOLAR Z_{out}) Jumper plug, W2, provides $Z_{out} \leq 1 \Omega$ or $\sim 93 \Omega$ for the rear-panel Bipolar output. Shipped in the $93-\Omega$ position.

BUSY/BUSY Jumper plug, W3, allows the Busy output to be a positive true or negative true logic signal. Shipped in BUSY (positive true) position.

PUR/PUR Jumper plug, W5, allows the Pile-Up Reject (PUR) output to be a positive true or negative true logic signal. Shipped in PUR (positive true) position.

INH/INH Jumper plug, W6, allows the INH IN input to accept either positive true or negative true logic signals. Shipped in INH (positive true) position.

TRI/GAUSS Jumper plug, W7, allows optimal livetime correction when used with ORTEC analyzers like the ADCAM® by connecting the BUSY output to the analyzer Busy In as described in Section 3.8. The jumper should be set to match the Unipolar Mode, TRI for Triangle and GAUSS for Gaussian. Shipped in TRI position.

4.6. DIFFERENTIAL INPUT MODE

When long connecting cables are used between the detector and preamplifier input, noise induced in the cable by the environment can be a problem. The differential input mode can be used with paired

cables from the preamplifier to suppress the induced noise.

BAL (DIFFERENTIAL INPUT GAIN BALANCE)

The BAL potentiometer is used to adjust the gain balance between the positive and negative inputs and to adjust the balance between the front- and rear-panel inputs when the differential (DIFF) input mode is used. The initial adjustment of Gain Balance is made by providing the same input to both the front- and rear-panel inputs. This can be accomplished by using a BNC "T" connector to feed the input signal on the front-panel input to the rear-panel input. Set the amplifier gain to maximum. Connect an oscilloscope to the unipolar output. While observing the signal on the oscilloscope, use a small screwdriver to adjust the Gain Balance (internal adjustment has been factory set, Fig. 4.5) potentiometer until the display on the oscilloscope shows minimum signal. Remove the BNC "T" connector when the adjustment is complete, and the positive and negative gains will be matched for use with NORM input.

If the differential input mode is being used, connect the differential input cable to the BNC connector on the rear panel. Adjust BAL potentiometer until there is minimum noise around the baseline of the output signal. If there is a problem in getting minimum noise, repeat the initial procedure with the BNC "T" and the adjustment.

DIFFERENTIAL INPUT SIGNAL The differential input signal or phantom is used only in the differential (DIFF) input mode. The normal preamp output is connected to the front-panel input with the amplifier input polarity set to match this signal. A second output cable must be added to the preamplifier with its center, signal pin connected to the preamplifier ground with the same value as the normal preamp output series resistor (usually 93.1Ω or 51Ω).

Many ORTEC preamplifiers have two Energy outputs, each with a $93.1-\Omega$ series resistor. For differential operation, one output is connected to the amplifier front-panel input. The second output is modified by connecting the preamplifier end of the series $93.1-\Omega$ resistor to ground within the preamp (soldering may be necessary). This second output

should be properly marked and connected to the rear-panel input. Both cables should be the same length and be run next to each other.

4.7. SYSTEM THROUGHPUT

To achieve the desired results in high-rate energy spectroscopy, the experimenter must consider not only the input rate, but also the unpiled-up output rate. The unpiled-up output rate is determined by the processing time of the shaping amplifier, the pile-up inspection time, and the input rate. For semi-Gaussian time-invariant filter amplifiers, the unpiled-up output rate is theoretically given by ²

$$r_o = r_i \exp(-T_D r_i) \quad (1)$$

where r_o is the unpiled-up output count rate, r_i is the input count rate, and T_D is the deadtime or effective processing time of the amplifier. The value of T_D is equal to the sum of the effective amplifier pulse width, T_w , and the time-to-peak of the amplifier output pulse, T_p . The type of deadtime in the shaping amplifier is referred to as extending deadtime since a second event arriving before the end of the initial deadtime ext ends the deadtime by an additional amplifier output pulse width, T_w from the occurrence of the second pulse.

A normalized plot of Equation (1) is shown as the solid line in Fig. 4.6. The maximum mean output rate equals $1/T_D \exp(1)$ and occurs when the mean input rate equals $1/T_D$. At this maximum output rate the deadtime losses are 63.2%. For input count rates exceeding $1/T_D$ the unpiled-up output rate decreases. When using a pile-up inspection circuit, the value of T_D is given either by the sum of T_w and T_p , or by the sum of T_p and the pile-up inspection time, whichever is larger.

Spectroscopy systems also have a deadtime that is caused by the digitizing time of the Analog-to-Digital Converter (ADC). This deadtime is a non-extending deadtime since events arriving during the

digitizing time are ignored. For non-extending deadtime the output rate is given by ²

$$r_o = \frac{r_i}{1 + r_i T_D} \quad (2)$$

where T_D is the digitizing time for the ADC and is designated T_M in Equation (3). This relationship is shown as the dashed line in Fig. 4.6. The maximum obtainable output count rate is $1/T_D$ and occurs at $r_i = \infty$.

When the ADC is considered as part of the spectroscopy system, the deadtimes of the amplifier and ADC are in series. The combination of the extending deadtime of the amplifier followed by the non-extending deadtime of the ADC is given by ²

$$r_o = \frac{r_i}{\exp[r_i(T_w+T_p)] + r_i[T_M - (T_w+T_p)] U[T_M - (T_w+T_p)]} \quad (3)$$

where $U[T_M - (T_w+T_p)]$ is a unit step function that changes value from 0 to 1 when T_M is greater than (T_w+T_p) . Equation (3) reduces to Equation (1) when T_M is less than (T_w+T_p) .

A plot of the unpiled-up amplifier output rate as a function of input rate for six values of shaping time is shown in Fig. 4.7. The measured deadtime, T_D is shown for each shaping time constant. The maximum value of the unpiled-up output rate increases with decreasing values of shaping time

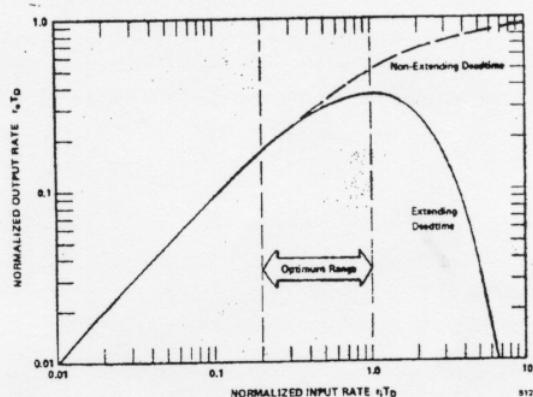


Fig. 4.6. Plot of Normalized Output Rate as a Function of Normalized Input Rate for Spectrometers with Simple Deadtime.

²R. Jenkins, R.L. Gould, and D.A. Gedcke, *Quantitative X-Ray Spectroscopy*, Marcel and Dekker, Inc., New York, (1980).

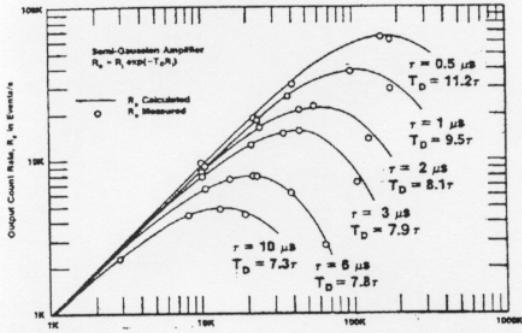


Fig. 4.7. Plot of the Unpiled-Up Amplifier Output Rate as a Function of Input Rate for Six Values of Shaping Time Constants.

constant. A set of throughput curves will remain nearly unchanged for a given amplifier for various energy ranges, detector types, and sizes.

The advantage of shorter shaping time constants to achieve higher output count rates is clearly shown in Fig. 4.7. However, shorter time constants also result in increased noise and increased charge collection time effects. Under worst case conditions, the noise increases inversely as the square root of the ratio of shaping time constants. The increase in the total energy resolution is the noise contribution combined in quadrature with the statistical contribution of the detector at the energy of interest. Consequently, the percentage of degradation in energy resolution can be much less than the percentage increase in noise.

4.8. CHARGE COLLECTION OR BALLISTIC DEFICIT EFFECTS

Charge collection distances in large-volume HPGe detectors are often 3 cm or more, resulting in charge collection times exceeding 300 ns.^{3,4,5} These charge collection times are due to the transit time of the holes and the electrons in germanium

³E. Sakai, "Charge Collection in Coaxial Ge(Li) Detectors," *IEEE Trans. Nucl. Sci.*, NS-15, 310, (1968).

⁴E. Sakai, T.A. McMath, and R.G. Franks, "Further Charge Collection Studies in Coaxial Ge(Li) Detectors," *IEEE Trans. Nucl. Sci.*, NS-16, 68, (1968).

⁵T.H. Becker, E.E. Gross, and R.C. Trammell, "Characteristics of High-Rate Energy Spectroscopy Systems with Time-invariant Filters," *IEEE Trans. Nucl. Sci.*, NS-28, 1, (1981).

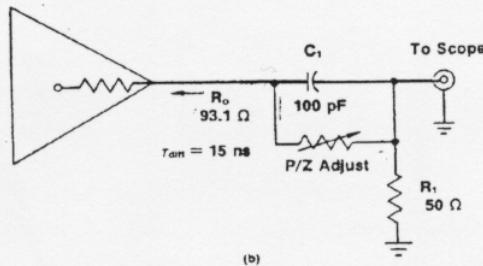
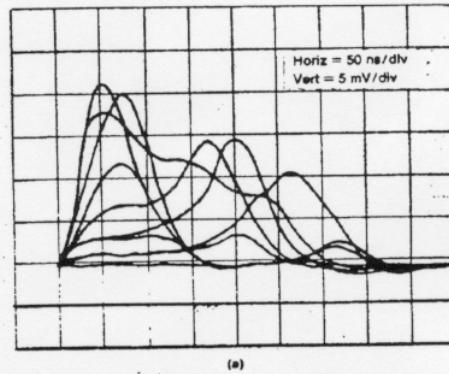


Fig. 4.8. Charge Collection Effect Waveforms. (a) Typical current Pulse Waveforms for a 28% Efficient HPGe Detector, and (b) the Simple Differentiation Circuit Used to Obtain the Current Waveforms.

and are not due to defects in the detector. Fig. 4.8(a) shows some typical current pulse waveforms from a 140-cm³ 28% efficient HPGe detector. These current pulse waveforms were obtained using the simple differentiation circuit shown in Fig. 4.8(b), which has a 15-ns time constant. The current pulses range in duration from 100 ns to greater than 350 ns. Pulses having equivalent total charge but different durations produce different output pulse heights when processed by a charge-sensitive preamplifier and a semi-Gaussian filter amplifier. This results in the distortion of the spectrum in direct proportion to the pulse amplitude or energy. This distortion is most pronounced at short shaping time constants. Figure 4.9(a) shows a portion of a spectrum obtained with a 2-μs shaping time, using the 1.33-MeV line of ⁶⁰Co. An equivalent spectrum using a 0.5-μs shaping time is shown in Fig. 4.9(b) and is significantly distorted.

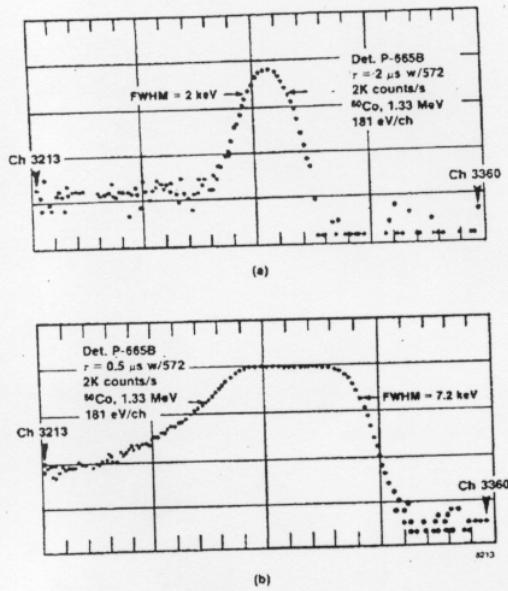


Fig. 4.9. Charge Collection Effect Spectrum. Logarithmic Display of Spectrum Taken with a 10% Efficient HPGe Detector for the 1.33 MeV ^{60}Co Line. (a) A 2- μs Shaping Time Constant and (b) a 0.5- μs Shaping Time Constant.

Charge collection time effects are of significant importance when using large-volume Ge detectors at high energy. The performance of two HPGe detectors is compared in Fig. 4.10 at two different energies. When using the 122-keV line of ^{57}Co , the principal cause of resolution degradation with decreased shaping time constant is the increase in noise. However, when using the 1.33-MeV line of ^{60}Co , the significant degradation in resolution is due to charge collection effects. The calculated resolution for the 10% detector at 1.33 MeV is shown as the dashed line in Fig. 4.10, and indicates approximately 2.0 keV FWHM at a 0.5- μs shaping time constant. The measured resolution under these test conditions was 7.2 keV, indicating that charge collection effects dominate. In Fig. 4.10, charge collection effects begin to appear at time constants less than 3 μs .

4.9. PILE-UP REJECTOR (PUR) AND LIVETIME CORRECTOR

An efficient pile-up rejector is incorporated in the amplifier to suppress the spectral distortion which is

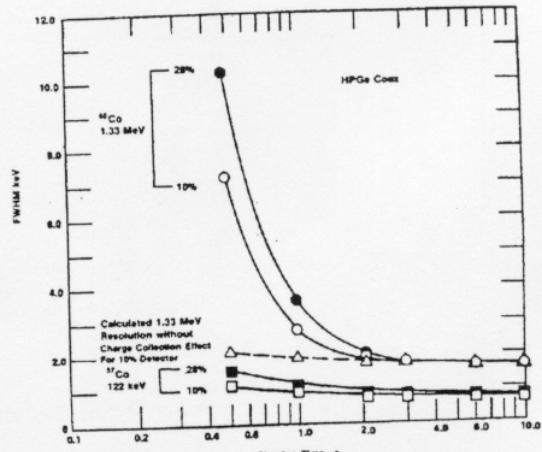


Fig. 4.10. Energy Resolution FWHM as a Function of Amplifier Shaping Time Constant for a 10% HPGe Detector and a 28% HPGe Detector for the 122-keV ^{57}Co Line and the 1.33-MeV ^{60}Co Line.

caused by pulses piling up on each other at high counting rates. High counting rate for pile-up is dependent on the dead time per pulse, T_B and hence the selected shaping time. T_B is 9 times the front-panel shaping time, T_c . High count rate for the PUR is when the normalized count rate $R T_i > 0.5$, where R is the amplifier input rate (see Fig. 4.6). For example, for 6- μs shaping R is 9 kHz and for 2- μs shaping, R is 28 kHz. Amplifier throughput for this condition using Equation (1) in Section 4.7 is 60% of the input rate. A multicolor pile-up rejector LED is included on the front panel to indicate the throughput efficiency of the amplifier. At low counting rates (pulse pile-up losses <40%) the LED flashes with a green color. At moderate counting rates the color changes to yellow. The color changes to red at high counting rates when the pulse pile-up losses are >70%.

The fast amplifier in the pile-up rejector includes a gated baseline restorer with its own automatic noise discriminator to eliminate the need for any operator adjustments. This function is also protected against negative overloads from pulsed reset preamplifiers. The PUR (pile-up reject) output logic pulse can be used at the gate or reject input of a multichannel analyzer to suppress pile-up in the recorded spectrum.

The block diagram for a gamma-ray spectroscopy system with pile-up rejection and live time correction is shown in Fig. 4.11.

FOR A RESISTIVE FEEDBACK PREAMP, CONNECT:

- Inhibit pulse from PUR to ADC PUR or ADC anticoincidence input.

- Livetime correction signal (Busy output) to the ADC Busy In.

ADDITIONAL CONNECTION FOR TRP (Transistor Reset Preamplifiers) Shown in dotted lines.

- Inhibit Output from TRP to the amplifier inhibit In.

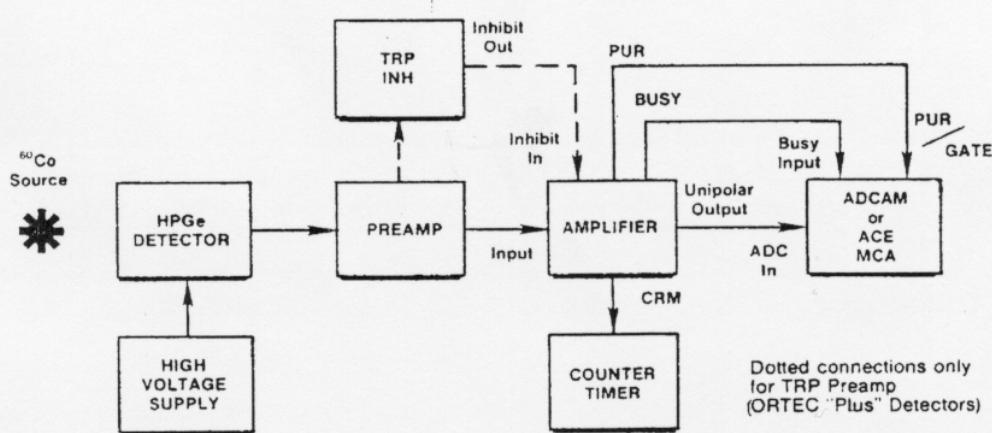


Fig. 4.11. Block Diagram for a Gamma-Ray Spectroscopy System with Pile-Up Rejection and Livetime Correction.

4.10. OPERATION WITH SEMICONDUCTOR DETECTORS

CALIBRATION OF TEST PULSER An ORTEC 419 Precision Pulse Generator, or equivalent, is easily calibrated so that the maximum pulse height dial reading (1000 divisions) is equivalent to a 10-MeV loss in a silicon radiation detector. The procedure is as follows:

- Connect the detector to be used to the spectrometer system, that is, preamplifier, main amplifier, and biased amplifier.
- Allow excitation from a source of known energy (for example, alpha particles) to fall on the detector.

c. Adjust the amplifier gain and the bias level of the biased amplifier to give a suitable output pulse.

d. Set the pulser Pulse Height control at the energy of the alpha particles striking the detector (e.g., set the dial at 547 divisions for a 5.47-MeV alpha particle energy).

e. Turn on the pulser and use its Normalize control and attenuators to set the output due to the pulser for the same pulse height as the pulse obtained in step c. Lock the Normalize control and do not move it again until recalibration is required.

The pulser is now calibrated; the Pulse Height dial reads directly in MeV if the number of dial divisions is divided by 100.

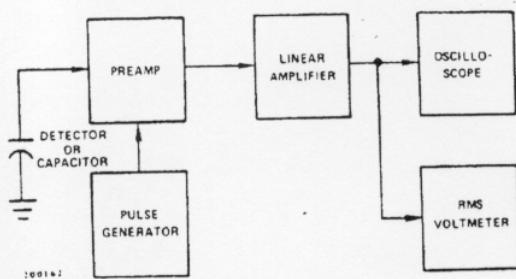


Fig. 4.12. System for Measuring Amplifier and Detector Noise Resolution.

AMPLIFIER NOISE AND RESOLUTION MEASUREMENTS As shown in Fig. 4.12, a preamplifier, amplifier, pulse generator, oscilloscope, and wide-band rms voltmeter such as the Hewlett-Packard 3400A are required for this measurement. Connect a suitable capacitor to the input to simulate the detector capacitance desired. To obtain the resolution spread due to amplifier noise:

- Measure the rms noise voltage (E_{rms}) at the amplifier output.
- Turn on the 419 precision pulse generator and adjust the pulser output to any convenient readable voltage, E_o , as determined by the oscilloscope.

The full-width-at-half-maximum (FWHM) resolution spread due to amplifier noise is then

$$N(FWHM) = \frac{2.35 E_{rms} E_{dial}}{E_o}$$

where E_{dial} is the pulser dial reading in MeV and 2.35 is factor for rms to FWHM. For average-responding voltmeters such as the Hewlett-Packard 400D, the measured noise must be multiplied by 1.13 to calculate the rms noise.

The resolution spread will depend on the total input capacitance, since the capacitance degrades the signal-to-noise ratio much faster than the noise.

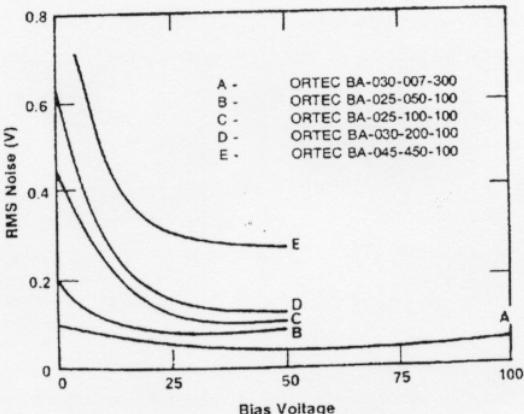


Fig. 4.13. Noise as a Function of Bias Voltage.

DETECTOR NOISE-RESOLUTION MEASUREMENTS The measurement just described can be made with a biased detector instead of the external capacitor that would be used to simulate detector capacitance. The resolution spread will be larger because the detector contributes both noise and capacitance to the input. The detector noise-resolution spread can be isolated from the amplifier noise spread if the detector capacity is known, since

$$(N_{det})^2 + (N_{amp})^2 = (N_{total})^2,$$

where N_{total} is the total resolution spread and N_{amp} is the amplifier resolution spread when the detector is replaced by its equivalent capacitance.

The detector noise tends to increase with bias voltage, but the detector capacitance decreases, thus reducing the resolution spread. The overall resolution spread will depend upon which effect is dominant. Figure 4.13 shows curves of typical noise-resolution spread versus bias voltage, using data from several ORTEC silicon surface-barrier semiconductor radiation detectors.

AMPLIFIER NOISE-RESOLUTION MEASUREMENTS USING MCA Probably the most convenient method of making resolution measurements is with a pulse height analyzer as shown by the setup illustrated in Fig. 4.14.

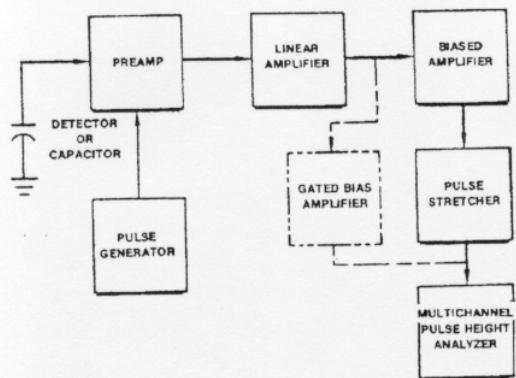


Fig. 4.14. System for Measuring Resolution with a Pulse Height Analyzer.

The amplifier noise-resolution spread can be measured directly with a pulse height analyzer and the mercury pulser as follows:

- Select the energy of interest with an ORTEC 419 Precision Pulse Generator. Set the amplifier and biased amplifier gain and bias level controls so that the energy is in a convenient channel of the analyzer.
- Calibrate the analyzer in keV per channel, using the pulser; full scale on the pulser dial is 10 MeV when calibrated as described above.
- Obtain the amplifier noise-resolution spread by measuring the FWHM of the pulser peak in the spectrum.

The detector noise-resolution spread for a given detector bias can be determined in the same manner by connecting a detector to the preamplifier input. The amplifier noise resolution spread must be subtracted as described in Section 4.10, "Detector Noise-Resolution Measurements." The detector noise will vary with detector size and bias conditions and possibly with ambient conditions.

CURRENT-VOLTAGE MEASUREMENTS FOR Si AND Ge DETECTORS The amplifier system is not directly involved in semiconductor detector current-voltage measurements, but the amplifier serves to permit noise monitoring during the setup. The detector noise measurement is a more sensitive method of determining the maximum detector

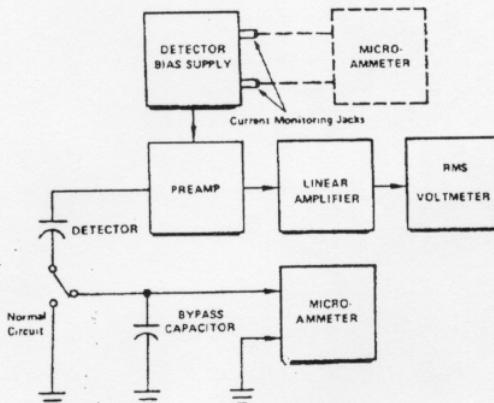


Fig. 4.15. System for Detector Current and Voltage Measurements.

voltage than a current measurement and should be used because the noise increases more rapidly than the reverse current at the onset of detector breakdown. Make this measurement in the absence of a source.

Figure 4.15 shows the setup required for current-voltage measurements. An ORTEC 428 Bias Supply is used as the voltage source. Bias voltage should be applied slowly and reduced when noise increases rapidly as a function of applied bias. Figure 4.16 shows several typical current-voltage curves for ORTEC silicon surface-barrier detectors.

When it is possible to float the microammeter at the detector bias voltage, the method of detector current measurement shown by the dashed lines in Fig. 4.15 is preferable. The detector is grounded as in normal operation, and the microammeter is connected to the current monitoring jack on the 428 detector bias supply.

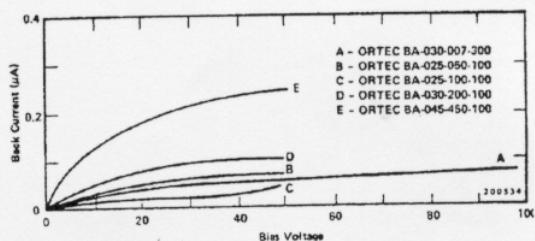


Fig. 4.16. Silicon Detector Back Current vs Bias Voltage.

4.11. OPERATION IN SPECTROSCOPY SYSTEMS

HIGH-RESOLUTION ALPHA-PARTICLE SPECTROSCOPY SYSTEM The block diagram of a high-resolution spectroscopy system for measuring natural alpha particle radiation is shown in Fig. 4.17. Since natural alpha radiation occurs only above several MeV, an ORTEC 444 Biased Amplifier is used to suppress the unused portion of the spectrum; the same result can be obtained by using digital suppression on the MCA in many cases. Alpha-particle resolution is obtained in the following manner:

- Use appropriate amplifier gain and minimum biased amplifier gain and bias level. Accumulate the alpha peak in the MCA.

b. Slowly increase the bias level and biased amplifier gain until the alpha peak is spread over 5 to 10 channels and the minimum- to maximum-energy range desired corresponds to the first and last channels of the MCA.

c. Calibrate the analyzer in keV per channel using the pulser and the known energy of the alpha peak (see Section 4.10, "Calibration of Test Pulser") or two known energy alpha peaks.

d. Calculate the resolution by measuring the number of channels at the FWHM level in the peak and converting this to keV.

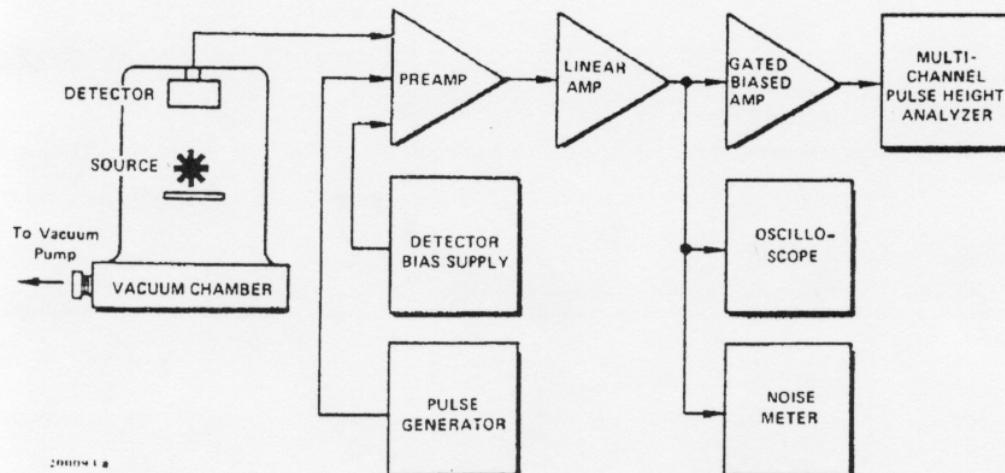


Fig. 4.17. System for High-Resolution Alpha-Particle Spectroscopy.

HIGH-RESOLUTION GAMMA-RAY SPECTROSCOPY SYSTEM A high-resolution gamma-ray spectroscopy system block diagram is shown in Fig. 4.18. Although a biased amplifier is not shown (an analyzer with more channels being preferred), it can be used if the only analyzer available has fewer channels and only higher energies are of interest.

When germanium detectors nitrogen cryostat are used, it is from about 1 keV FWHM up that are cooled by a liquid possible to obtain resolutions to 4 keV (depending on the energy of the incident radiation and the size and quality of the detector). Reasonable care is required to obtain such results. Some guidelines for obtaining optimum resolution are:

- a. Keep interconnection capacities between the detector and preamplifier to an absolute minimum (no long cables).
- b. Keep humidity low near the detector-preamplifier junction.
- c. Operate the amplifier with the shaping time that provides the best signal-to-noise ratio.
- d. Operate at the highest allowable detector bias to keep the input capacity low.

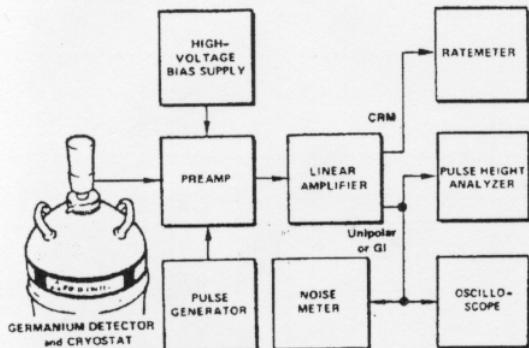


Fig. 4.18. System for High-Resolution Gamma Spectroscopy.

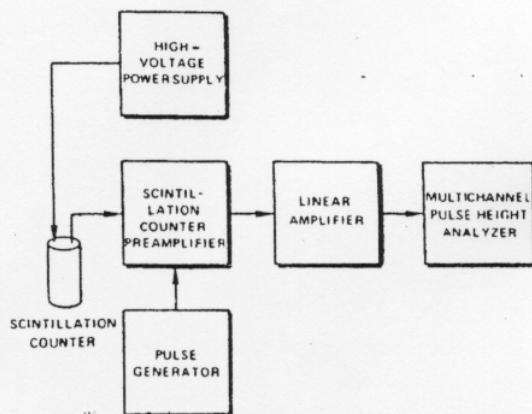


Fig. 4.19. Scintillation-Counter Gamma Spectroscopy System.

SCINTILLATION-COUNTER GAMMA SPECTROSCOPY SYSTEMS The ORTEC 671 can be used in scintillation-counter spectroscopy systems as shown in Fig. 4.19. The amplifier shaping time constants should be selected in the region of 0.5 to 1 μ s for NaI or plastic scintillators. For scintillators having longer decay times, longer time constants should be selected.

X-RAY SPECTROSCOPY USING PROPORTIONAL COUNTERS Space charge effects in proportional counters, operated at high gas amplification, tend to degrade the resolution capabilities drastically at x-ray energies, even at relatively low counting rates. By using a high-gain low-noise amplifying system and lower gas amplification, these effects can be reduced and a considerable improvement in resolution can be obtained. The block diagram in Fig. 4.20 shows a system of this type. Analysis can be accomplished by simultaneous acquisition of all data on a multichannel analyzer or counting a region of interest in a single-channel analyzer window with a counter and timer or counting ratemeter.

4.12. OTHER EXPERIMENTS

Block diagrams illustrating how the 671 and other ORTEC modules can be used for experimental setups for various other applications are shown in Figs. 4.21 through 4.24.